

Product Name: Halo Silicone Breastmilk Collector

Materials: Collectors - Medical-Grade Silicone
Tube Cleaning Brush - Nylon, Stainless Steel

Dimensions: 1/2 x Breastmilk Collector(s) - 103 x 101 x 56mm
1 x Tube Cleaning Brush - 8 x 8 x 99mm

Colour: Clear

Capacity: 3 oz./90ml

Product Description:

Soft suction breast pump, milk collector - why not combine both with our Halo Silicone Breastmilk Collector! This nifty little device acts as a manual wearable breast pump, making it ideal for those who prefer lighter suction on days you need to pump but have things to do. Simply press the back and attach it to your breast, pop on your bra, and get on with your day! Alternatively, if you just want to save your leakage, simply put it into your bra without any suction and passively collect. It's a win-win for you, no matter your needs. The Halo stands stably between uses to prevent any milk spills, and the curved easy-pour opening lets you transfer your precious liquid gold into a bottle or storage container without any mess. The easy-to-read measurements ensure you can track what you've collected at a glance. Easy to use and clean - it even comes with a brush to make washing a breeze. Halo is the perfect choice for mums who need a light-suction pump and collector in one!

Features:

- Doubles as a wearable light-suction breast pump and breast pad alternative.
- Simply press the back to activate gentle continuous natural suction for easy milk expression or collection.
- Discreet, hands-free design ensures comfortable use anytime, anywhere.
- Curved easy-pour opening for spill-free milk transfer.
- Clear measurement markings help easily track how much milk you've collected.
- Guide groove with indication arrow clearly indicates milk transfer.
- Concave shape protects your nipples and naturally moulds to your breast for a secure, comfortable fit.
- Stable base allows the collector to stand upright, preventing spills after use.
- Made from premium 100% medical-grade silicone that is safe for your skin and breastmilk.
- Easy-care, one-piece design with dedicated cleaning brush - simply boil to sterilise.
- Portable and travel-friendly - no cords or batteries required.
- BPA, PVC, and phthalate-free.


Cleaning and Care:

Wash before and after each use. Although dishwasher safe, we recommend hand washing with warm, soapy water and rinsing thoroughly. To clean, use the included cleaning brush to scrub thoroughly with warm, soapy water and rinse well. To sterilise, use a steam steriliser or boil in water for 2-3 minutes. Brush NOT suitable for microwaves, sterilisers or boiling water due to stainless steel and nylon components. Only soft bristle brushes or soft sponges should be used to clean this product as hard scourers may scratch the surface. Do not use any bleach-based cleaners or tablets to clean or sterilise your Haakaa products. Store product in a cool, dry place and avoid direct sunlight.


Caution:

- Check product condition regularly. Replace if this product shows any signs of damage.
- Adult use only.
- This product is not a toy. Do not use this product for anything other than its intended use.
- Keep out of reach of children and pets when not in use.
- Do not store near any sharp objects.
- Your silicone products may grow cloudy after repeated use, especially in hard water areas. This is normal and will not affect the performance of your product.
- Due to stainless steel and nylon components, do not place tube cleaning brush in microwaves, sterilisers or boiling water.


Instructions:

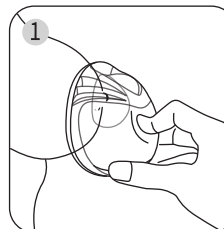

Press down on the back of the collector and centre the main opening over your nipple.

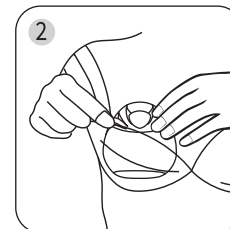


Press to your breast and release to create suction. Secure in your bra.

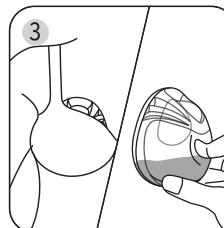

To restart suction, take out your collector and press the back again, taking care not to spill any collected milk.

You may need to transfer your milk to a clean, sterilised container before restarting.

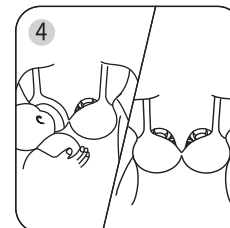

To use with no suction, pop into your bra without pressing the back. Wear one on the other side to save let-down while nursing, or wear two at a time to express or collect on-the-go.

PLEASE NOTE: Volume is an approximation and should be used for general reference only. Empty milk collector regularly and store in accordance with your country's guidelines on the safe storage of breastmilk. Milk collected over extended periods and left at room or body temperature for too long should not be fed to your baby. The Halo utilises lighter suction strength. If you require stronger suction, please check Haakaa's range of manual silicone breast pumps.