

Product Name: 360° Maternity Lift & Support Band

Materials: Band - 49% Nylon, 33% Polyester, 18% Rubber (Elastodiene)
Support Stays - PP

Dimensions: M - 130 × 20cm
L - 150 × 20cm

Colour: Grey

Sizes: M/L

Product Description:

Ease the strain of pregnancy with our 360° Maternity Lift & Support Band, designed to cradle your growing belly while relieving lower back pressure. This one-piece ergonomic band features upper and lower straps that form a secure 360° wraparound, offering balanced support to maintain and improve your posture. The wide panel provides enhanced waist support, distributing weight more evenly and reducing discomfort during extended periods of standing or walking. Three removable stays let you choose between firm reinforcement and relaxed flexibility for amazingly comfortable customisation. High-quality hook & loop tape ensures easy wear and removal. Made from soft, breathable fabric, this support band provides a snug, comfortable fit that moves with you throughout the day. It's the ideal choice for expectant mothers seeking reliable support in the latter stages of pregnancy!

Features:


- O-shaped wraparound provides 360° belly and back support.
- Wide panel stabilises your spine and distributes belly weight evenly.
- Effectively relieves lower back, pelvic, and hip pain during pregnancy.
- Removable stays strengthen support when inserted, and offer a more flexible feel when removed.
- One-piece design enables easy adjustment, convenient for late pregnancy.
- Breathable bamboo charcoal fabric ensures lightweight, all-day comfort.
- Ultra-soft, smooth binding for flexible use without chafing.
- Contoured hip design prevents shifting, allowing you to move freely or sit comfortably without pressure.
- Fully adjustable straps provide personalised support from the second to third trimester.
- Targeted support keeps pregnant bodies stable during daily activities, at home or on-the-go.
- Designed with two size options to provide optimal support for a range of body types.
- Stays durable and elastic after repeated washing, providing lasting support.

Cleaning and Care:


Remove support stays and clean this product regularly. Hand wash in warm water (under 30°C) and rinse thoroughly. Lie flat to dry. Do not forcefully scrub this product. Hard scourers may scratch the surface. Not suitable for chlorine, bleach, ironing, or dry cleaning. Do not machine wash or tumble dry. Methods of cleaning or use not listed in this manual may damage product. Store in a cool, dry place and avoid direct sunlight when not in use.

PLEASE NOTE: Always consult your physician before using this product, as wear times will vary depending on your individual situation. We highly recommend wearing the band over your clothes to minimise the risk of chafing or discomfort. Band is not intended for waist training purposes.


Instructions:

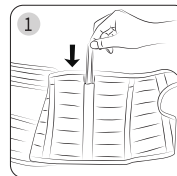


Insert stays if desired. Place the longest stay in the centre slot and the two shorter stays in the slots on either side. Insert with the raised side of the stays facing outward.

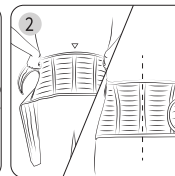

Align the centre of the back of the band with your spine, ensuring the widest part is positioned directly along your spine.

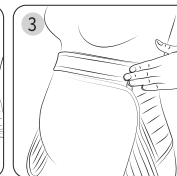

Wrap the upper straps around to the top of the abdomen, securing them to provide support and coverage for the upper belly.

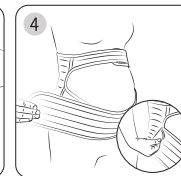

Wrap the lower straps around the bottom of the abdomen and secure at the front. To remove, simply pull to release.

NOTE: Band may shift slightly while wearing - reposition as needed. Ensure the stays are inserted securely into the correct slots.

BAND SIZING GUIDE

	M	L
Recommended Upper Belly Circumference	≤94cm / ≤37 inches	≤98cm / ≤38.6 inches
Recommended Lower Belly Circumference	≤118cm / ≤46.5 inches	≤140cm / ≤55.1 inches

NOTE: Please measure your lower belly below your navel, and your upper belly just below your underbust.

Caution:

- We do not recommend placing this product directly onto your skin. If you choose to do so, avoid wearing it for long periods of time to minimise the risk of chafing or discomfort.
- Always consult your medical professional prior to use. This product is not suitable for mums with minimal belly growth. Please assess your own body condition to determine the appropriate time to start using it.
- Do not wear if you have an allergy to any of the materials listed.
- Adult use only. Keep out of reach of children and pets.
- Check product condition regularly. Replace if product shows signs of damage.
- This product is a pregnancy assistive tool and is not a medical device.
- We do not recommend wearing this product for extended periods of time or while sleeping.
- Discontinue use if you feel any discomfort.
- Consult your healthcare professional if you're unsure about the correct fit during pregnancy. Ensure proper wear and avoid wearing too tightly for optimal comfort and support. Not to be used during labour.
- Avoid vigorous activities or excessive bending while wearing this product to prevent discomfort.