尺寸:
210x225mm

工艺要求:
70g双胶纸

半成品号:
E10949

haakaa® LATCHFLOW SYSTEM™ - 20 mL DISPENSER Product Manual English

BHK185 E10949

Product Name: LatchFlow System™ - 20 mL Dispenser

Materials:
Nipple Shield - Medical-Grade Silicone
Feeding Tubes, Dispenser - Food-Grade Silicone
Extension Tubes - Food-Grade Silicone, PBT
Tube Connector, Storage Frame & Case - PP
Cleaning Brush - Stainless Steel, Nylon, PP
Medical Tape - Nonwoven Fabric

Dimensions:
Nipple Shield (S) - 131 x 63 x 77 mm
Nipple Shield (M) - 131 x 63 x 80 mm
Nipple Shield (L) - 131 x 63 x 81 mm
Dispenser - 39 x 39 x 69 mm
Extension Tubes - 9 x 9 x 390 mm
Feeding Tubes (S) - 11 x 10 x 400 mm
feeding end diameter: 1.6 mm
Storage Frame - 79 x 69 x 49 mm
Cleaning Brush - 16 x 9 x 150 mm
Storage Case - 180 x 128 x 60 mm
Medical Tape - 9100 x 13 mm

Volume:
Dispenser Capacity: 0.7 fl oz (20 mL)
Connection Opening: 3.8 mm

Nipple Shield Standard: EN 14350

Intended Use:

This system is intended to provide supplemental milk to babies while they feed at the breast, using either the nipple shield or a feeding tube. Use only as recommended by a lactation professional.

Product Description:

Designed for at-breast supplementation, this dispenser-based feeding system helps provide additional milk while baby remains at the breast. The soft 20 mL dispenser allows gentle, adjustable milk flow with a simple squeeze, while the nipple shield and feeding tubes support flexible feeding options.

The direct-flow dispenser features hand-held squeeze control for precise, responsive flow. The easy-to-connect extension and direct-feed tubes offer flexible feeding options for at-breast supplementation, and the form-holding frame helps the shield maintain its shape for neat, hygienic storage. Each component has been carefully designed to make this system easy to assemble, use, clean, and sterilise.


Features:

- Provides extra supplemental milk to baby while supporting continued breastfeeding, with flexible direct or shield-assisted feeding options.
- Helps mothers continue breastfeeding and maintain skin-to-skin contact while lactation is being established.
- Perfectly-sized for smaller, more frequent feeds.
- Suitable for colostrum, expressed milk, or formula as needed.
- Double-layered silicone shield with built-in channels helps guide supplemental milk toward the nipple.
- Ultra-thin, medical-grade silicone shield has unique cut-outs to allow closer contact for a more natural feeding experience.
- Nipple shield fits Haakaa LatchFlow System™ dispensers and pouches, other supplemental feeding systems and oral feeding syringes via the 3.8 mm opening (only the dispenser is included).
- Openings on the nipple shield allow milk to pass through from mum's nipple.
- Easy-squeezy 20 mL dispenser features one-handed milk flow control.
- Extra tube for optional extension use or as a replacement to meet every need.
- Includes a shape-retaining frame for the nipple shield and a storage case to keep all parts neatly organised.
- Minimal parts to set up and disassemble – safe in dishwashers, sterilisers and boiling water.
- BPA and phthalate-free.

Cleaning and Care:

- Before first use, disassemble and clean all parts of this product before placing silicone components in boiling water for 5 minutes to sterilise. Boil plastic parts, including dispenser cap & tube connectors, for 2-3 minutes to sterilise (do not boil any longer as this may impact the lifespan of plastic parts). Wash before and after each use.
- To clean the tubing, fill dispenser with warm soapy water and connect the feeding tube as normal. Turn the dispenser upside down and squeeze to let water flow through the tubing. Squeeze several times to flush through. Squeeze and roll tube between your fingers to ensure thorough cleaning of residue.
- To rinse the tubing, fill dispenser with clean drinking-quality water and connect the tube. Turn dispenser upside down and squeeze the water through the tube. Repeat as necessary until thoroughly rinsed. If any water is trapped inside the tubing after cleaning, drain the water from dispenser, connect the tube then squeeze to flush it out. Clean tubing as soon as possible after feeding to prevent milk drying inside. It is important all residual moisture dries.
- To clean the shield, press the centre dome through from the other direction to separate it from the outer flange. Use the included cleaning brush to gently scrub all shield surfaces, both inside and out. Leave inverted to thoroughly dry.
- To clean other parts of the product, although dishwasher safe, we recommend hand washing in warm soapy water, rinsing and drying thoroughly.
- To sterilise all parts except the cleaning brush, boil in water for 2-3 minutes. Brush **NOT** suitable for microwaves, sterilisers or boiling water due to stainless steel components. Do not use any bleach-based agents or sterilising tablets to clean this product. Do not use a UV steriliser to clean this product as it may impact its lifespan. When boiling to sterilise, ensure parts do not come into direct contact with the sides or bottom of the container, as this may cause deformation.
- Only soft bristle brushes or soft sponges should be used to clean this product, as hard scourers may scratch the surface. Store in a cool, dry place and avoid direct sunlight.

Find Your Perfect Fit:

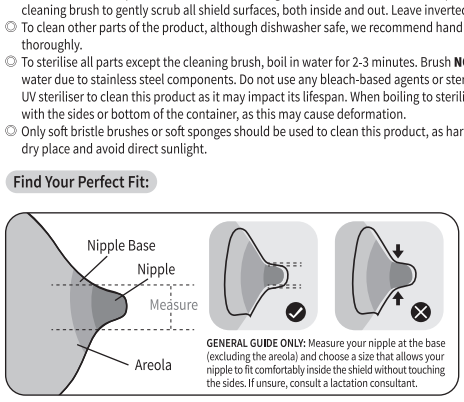


GENERAL GUIDE ONLY: Measure your nipple at the base (excluding the areola) and choose a size that allows your nipple to fit comfortably inside the shield without touching the sides. If unsure, consult a lactation consultant.

Shield Opening Size	Recommended Nipple Size
S: 17 mm	13-16 mm
M: 20 mm	16-19 mm
L: 24 mm	19-23 mm


Instructions:

Ensure each part is correctly fitted and securely connected before use.

Feeding from direct tube:

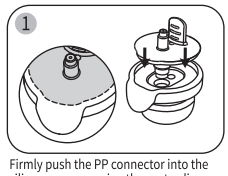


Firmly push the PP connector into the silicone cap, ensuring the protruding end of the connector is aligned with the matching opening in the cap. Check there are no gaps.

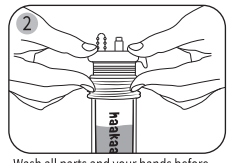


Wash all parts and your hands before use. Fill dispenser, pouch* or feeding syringe* with supplemental milk. Push the silicone cap and connector firmly into the pouch opening.

* Not included

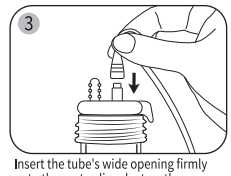


Insert the tube's wide opening firmly onto the protruding duct on the connector.

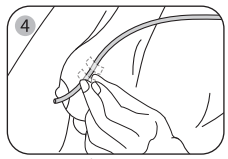

Place the end of the tube 0.5 cm away from your nipple and secure with the included medical tape.



Introduce your baby to your nipple. Place the end of the tube gently into their mouth to begin feeding.

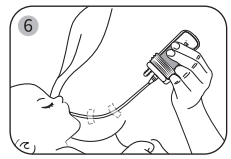


Gently squeeze the silicone dispenser to start and adjust the milk flow when operating.


WARNING: CHOKING HAZARD – SMALL PARTS. Keep out of reach of children under 3 years.

PLEASE NOTE: Tubes, connector and cap come assembled. In the event that reassembly is needed, please refer to the assembly instructions listed above. Always follow your country's guidelines and recommendations for safely storing breastmilk.
If using the nipple shield with other supplemental feeding systems or oral feeding syringes (not included), ensure the connector fits securely into the 3.8 mm opening before use.


Feeding from nipple shield:

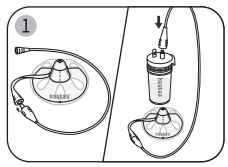


Fill the included dispenser, or a compatible pouch or oral feeding syringe (not included), with supplemental milk. Ensure the container and tubes are properly and securely assembled.

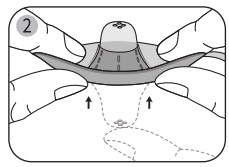


Ensure the two layers of the shield overlap correctly. Carefully turn brim upwards.

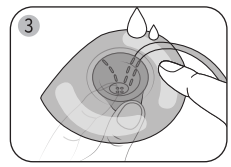


Place a few drops of drinking-quality water or breastmilk on the underside of the brim to help the shield stick to your skin.

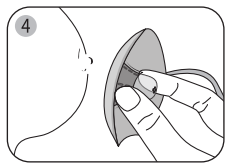

Centre the shield over your nipple, enclosing as much of your nipple as possible.

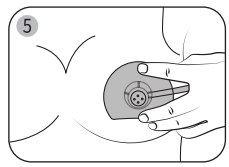

Press the shield down onto your skin to turn the brim back down. The base should sit flush against your skin.

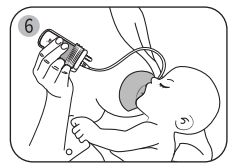

Gently squeeze the silicone dispenser to start and adjust the milk flow when operating. NEVER squeeze forcefully to avoid milk spurting out.

NOTE: If using the nipple shield with LatchFlow System™ pouches, other supplemental feeding systems or oral feeding syringes (none included), ensure the connector fits securely into the 3.8 mm opening before use.

Caution:

For your child's safety and health
WARNING!

- Continuous and prolonged sucking of fluids will cause tooth decay.
- Always check food temperature before feeding.
- Throw away at the first signs of damage or weakness.
- Keep components not in use out of the reach of children.
- Never attach to cords, ribbons, laces or loose parts of clothing. The child can be strangled.
- Always use this product with adult supervision. Only adults should attempt to feed babies from this product. Secure the feeding tube during use and ensure it does not wrap around the baby's neck.
- This is a single-user product. Use by more than one user without adequate cleaning may present a health risk and should be avoided.
- We recommend replacing nipple shield regularly.
- This product should only be used under the advice of a lactation professional. Always consult a medical professional for further advice.
- Check all tubing, connectors and parts are securely assembled before each use.
- Never use feeding teats as a soother.
- Do not leave a feeding teat in direct sunlight or heat, or leave in disinfectant ("sterilising solution") for longer than recommended, as this may weaken the teat.
- This product is not a toy. Do not use this product for anything other than its intended use.
- Do not warm breastmilk or infant formula in a microwave oven as this may cause uneven heating and could scald the infant.
- If using third-party feeding components, only use ones that securely fit the 3.8 mm opening. Do not force incompatible parts into the opening.
- Do not forcefully squeeze or shake at any time.
- The cap is designed with a hole to direct milk. It is not 100% leakproof for milk storage.
- Do not fill dispenser beyond top measurement marks.